


BioMicroLab

LabelPro

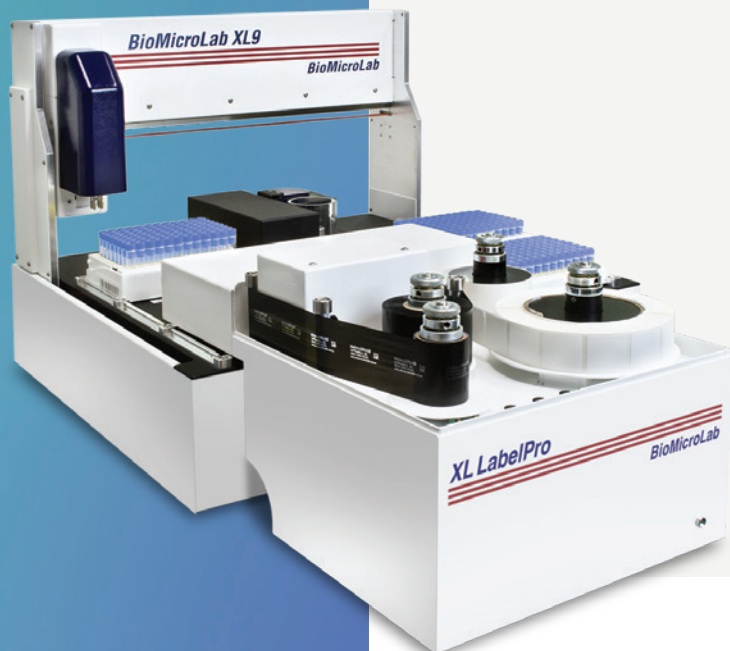

XL9 LabelPro automated labeling of microtubes

The XL9 tube handler with LabelPro automates tube labeling and can be used exclusively as a labeler, or integrated with other tasks in complex workflows.

automated tube labeling

- Robotics label, re-array and sort microtubes
- Label design software embedded within the tube handling software makes both designing labels and running labeling projects straight-forward.

benefits

- Easily affix human readable labels to labware
- Eliminate the headaches associated with manual tube labeling
- Fast start-up and unattended operation
- All electric operation; no compressed air
- Small footprint with convenient design

versatility for your laboratory

- The labeling module can be removed when not in use to free up 2 deck positions
- The labeler can run in Standalone Mode to label larger or various sized vials

label features

- Print 1D & 2D barcodes, logos, symbols and human readable text on white or any color label, using any font on your computer
- Print identical or unique labels
- 300 – 600 dpi thermal printer
- Label size typically 0.4" x 1.3" or 0.6" x 1.3"
- Labels designed for liquid nitrogen (-196°C) and freezer (-80°C) environments

compatible with 96-format microtubes

- Compatible with ANSI/SLAS 96-format microtubes from major manufacturers such as FluidX, Thermo Scientific, Matrix, Micronic, ABgene, Corning and more
- Vials up to 1" (25mm) diameter in standalone mode



**Color labels
now available!**

 **sptlabtech**

XL9 LabelPro

labels and printing ribbon

- Labels designed for liquid nitrogen (-196°C) and freezer (-80°C) environments
- Thermal printer ribbon resin is rated super premium for DMSO and alcohol chemical resistance (T84/T85)
- Wide variety of label material including plastic/permanent acrylic (CryoLabel®)
- 1,200 – 2,000 labels per roll (varies with label size)

speed

- Label up to 180 vials per hour [20-25 sec / tube]

Speeds vary with additional automation operations like volume detection, re-arranging, barcode decoding, etc.


software

Design and layout of labels is accomplished with our label design software. Labeling functions are embedded in the XL Work List Manager Software which provides the ability to print the same data on all labels or merge variable data on each label. Project based workflow allows for easily reprinting and applying commonly run labels.

- User-friendly Windows-based software
- Design your label with text and/or barcodes
- Input worklist file controls all tube operations (.csv)
- Merge data from external sources (.csv or .xml)
- Easy to repeat previous labeling jobs
- Output file is user-definable (.csv)
- LIMS integration ready



LabelPro™ Model Comparison

for microtubes		for microtubes & vials		
XL9 LabelPro™	XL20 LabelPro™	XL100 LabelPro™	XL200 LabelPro™	
sizes and configurations for every lab				
compatible labware & labels				
plastic or glass 96 format microtubes			✓	
Cryovials, 15mL conical tubes, 4mL, 8mL glass vials	--	--	✓ up to 125mm high	
Typical Label Sizes	0.4" x 1.3" label typical for Matrix/Micronic 1.4mL microtubes		0.4" x 1.3" label typical for Matrix/Micronic 1.4mL tubes 1"x1" labels typical for Cryovials, 2mL vials	
Print Human Readable Text and 1D & 2D Barcodes				
rack capacity by removing the LabelPro™ when not in use, 2 additional rack positions are available.				
SLAS/ANSIRack Capacity	up to 7 tube racks (672 microtubes)	up to 18 tube racks (1728 microtubes)	up to 18 tube racks (1728 microtubes; 864 tubes/vials)	up to 28 tube racks (2688 microtubes; 1344 tubes/vials)
SLAS/ANSI Rack Capacity with balance for volume detection	up to 5 tube racks (480 microtubes)	up to 16 tube racks (1536 microtubes)		

The LabelPro™ labels vials up to 1" (25mm) in diameter when not attached to a tube handler.